




THD BANDY

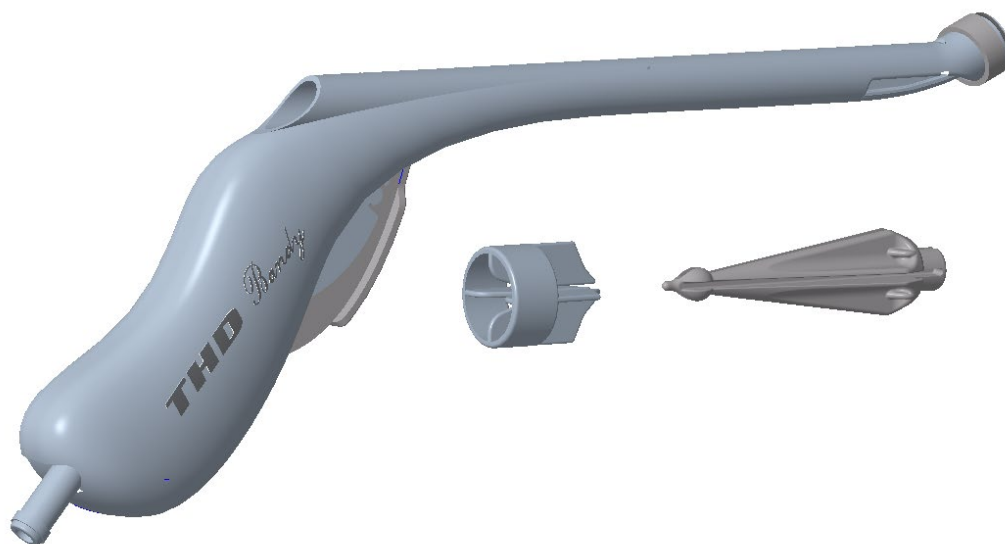

Use

THD Bandy is a disposable suction ligator for the treatment of hemorrhoidal disease.

General Description

THD Bandy is composed of an elastic suction ligator and a ring loading system, ready for use.

THD Bandy connects to a suction system, to control the quantity of aspirated tissue, around which the elastic band will be released by pulling the trigger.

The special design of the handle allows use with one hand as well as a fast and secure elastic release.

Features

- Made of a medical grade Polycarbonate
- Complies with the biocompatibility requirements of the ISO 10993 standard
- Supplied with disposable rings
- Loading with one hand using loading cone, included
- Ready to use
- Lightweight for easy handling
- FDA Regulatory Class II
- Individually packaged
- Latex Free



Advantages

- Can be used with one hand
- Easy loading of the elastic

Dimensional Features	
Length except handle (A)	150 mm (5.90 in)
Internal opening diameter (B)	10 mm (0.39 in)
Suction connection diameter (C)	8 mm (0.31 in)

THD Bandy			
REF	Q.ty	Description	
800183-25	25	THD Bandy	
<i>include</i>	/	25	Disposable suction ligator
<i>include</i>	/	25	Cones for ring loading
<i>include</i>	/	25	Ring loader
<i>include</i>	/	50	Rubber band ligator